

Solido 3

Universal bed server



Patient

- A tilting feature for reading
- A more pleasant stay on the bed



Staff

- Very easy table controllability
- Ergonomic table handling



Management

- A modern design with colour options
- Highly functional at an affordable price



Solido 3

The Solido 3 serving table is a practical accessory that will make the patient's stay on the bed more pleasant. It has a number of ergonomic features that makes handling very easy and intuitive for both the patient and the nursing staff.

Simple control

The Solido 3 server can be equipped with a spring for upward direction control without the use of the control handle.

Practical tilting

A tiltable eating board will be appreciated by patients when reading or using a tablet.

Ergonomic control

Thanks to the ergonomic features and special column structure, handling of the table is very easy.

- 1 A handle for easy handling.
- 2 Thanks to the increased space, the patient can pull the table even closer to the body which supports easier eating.
- 3 Exceptionally easy adjustment of the table height with a lever.
- 4 The integrated stopper protects books and other items from falling, even when the eating board is tilted.
- 5 The metal gallery serves as a practical stopper for items placed on the table.






HandyClip

The unique shape of the column allows you to attach a clip that can be shifted vertically. Various accessories can be attached to the clip – a water bottle holder, a urine container hook, a small shelf for a phone, etc. The patient always has everything he/she needs within arm's reach.



Design: functionality and aesthetics

The eating board has clean lines for easy cleanability, and it is produced in several colour schemes. The colour of the table can thus be matched to that of the bed and room amenities.

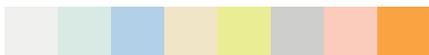

Useful details

- The smooth column, seamless technology and practical structural design make it easy to clean the T-shaped table.
- The T-shaped undercarriage can easily fit under the bed even in case of a fifth castor.
- Choose between:
 - The low clearance undercarriage with two small integrated castors.
 - The C shape undercarriage with the two larger 50 or 75 mm castors to move the table more easily and facilitate its maintenance.
- Table load-bearing capacity 30 kg.



Technical parameters

Inner dimensions of eating board	Fixed part	11×35 cm
	Big part	64×35 cm
External dimensions of eating board		94.5×45 cm
Height adjustment (depends on the castors and undercarriage)	C undercarriage and 50 mm castors	103×71 cm
	C undercarriage and 75 mm castors	104.5×72.5 cm
	T undercarriage	102×70 cm
Load-bearing capacity	In middle of eating board	30 kg
Dimensions of undercarriage	C undercarriage and 50 mm castors	80.5×44 cm
	T	79.5×43.5 cm
Height of undercarriage	C undercarriage and 50 mm castors	11.3 cm
	T	3.4 cm
Tilting angle		0° - 8° - 21° - 31° - 40° - 48°



Members of LINET Group

LINET spol. s r. o.

Želevčice 5 | 274 01 Slaný | Czech Republic

tel.: +420 312 576 400 | fax: +420 312 522 668 | e-mail: info@linet.com | www.linet.com

