

HELPING CONNECT BEDSIDE PRACTICE TO LEADERSHIP STRATEGY TO HELP REDUCE THE RISKS OF CLINICAL PRACTICE VARIATIONS

Being a caregiver is rewarding and challenging at the same time given the amount of multi-tasking that is required. We understand this and are here to assist the caregivers with their workflow by providing LINIS SafetyPort - a solution for automated export of data from LINET beds (e.g. safe position of the bed, bed exit alarm, patient weight).



NURSES

- 24/7 data overview of all the beds
- Data evaluation of patient mobilization
- Data to help reduce and manage pressure ulcers



MANAGEMENT

- Overview of bed features utilization
- Actionable data aggregated in graphs and charts
- Detailed overview of clinical practice at your facility



SERVICE & MAINTENANCE

- Real time data about the bed usage
- Complete history of error codes
- Bed localization in a hospital



PATIENT SAFETY

16 data points:

- controlling patient's weight and presence
- monitoring safety position of the bed



RESPIRATORY & MOBILIZATION

12 data points:

- improving mobilization practices
- preventing pressure injuries



BED MANAGEMENT

7 data points:

- monitoring location of bed and important maintenance data
- providing overview of its equipment



UTILIZATION

4 data points:

- providing hospital with actionable data overview of bed features

DASHBOARD

Customizable columns

History of all beds

Administration Interface

LOCATION	NAME	NOTE	BED MONITORING STATUS	OUT OF BED STATUS	SAFETY POSITION	SIDERAILS	BRAKES	LOWEST POSITION	BACKREST	TRENDELENBURG	LATERAL TILT	PATIENT WEIGHT
101 [A]	S. Z.	Disabled patient, paralyzed below...	ON - Inoner	IN	✓	✓	✓	✓	30°	0°	0°	58,7 kg
102 [B]	K. L.	Patient needs assistance	ON - Outer	OUT	✗	✗	✓	✓	0°	ATX 10°	0°	67 kg
103 [A]	J. M.	After surgery, latex allergy	OFF	IN	✗	✓	✓	✓	10°	0°	L 30° ALT	89 kg
104 [B]	A. N.	Patient rapidly losing weight	ON - Inoner	IN	✓	✓	✓	✓	30°	0°	0°	58,7 kg
105 [A]	R. Z.	Patient needs assistance	OFF	OUT	✗	✗	✓	✗	0°	0°	0°	67 kg
106 [A]	Not assigned		ON - Inoner	IN	✓	✗	✗	✗	30°	ATX 10°		
107 [B]	M. N.	Patient with diabetes	OFF	IN	✓	✓	✓	✓	30°	0°		
108 [A]	A. K.	Patient rapidly losing weight										Disconnected
109 [B]	Z. L.	Patient needs assistance										Disconnected
Unknown	R. Z.	Patient needs assistance	OFF	OUT	✗	✓	✓	✓	0°	0°		
Unknown	Z. L.	Patient needs assistance										Disconnected

24 hours View Care Notes Mobilization Statistics Insight

Detail of each bed

COMPATIBLE BEDS

22+
data points



Eleganza 3

22+
data points



Eleganza 4

37+
data points




Eleganza 5

39+
data points



Multicare, Multicare LE

CONNECTIVITY

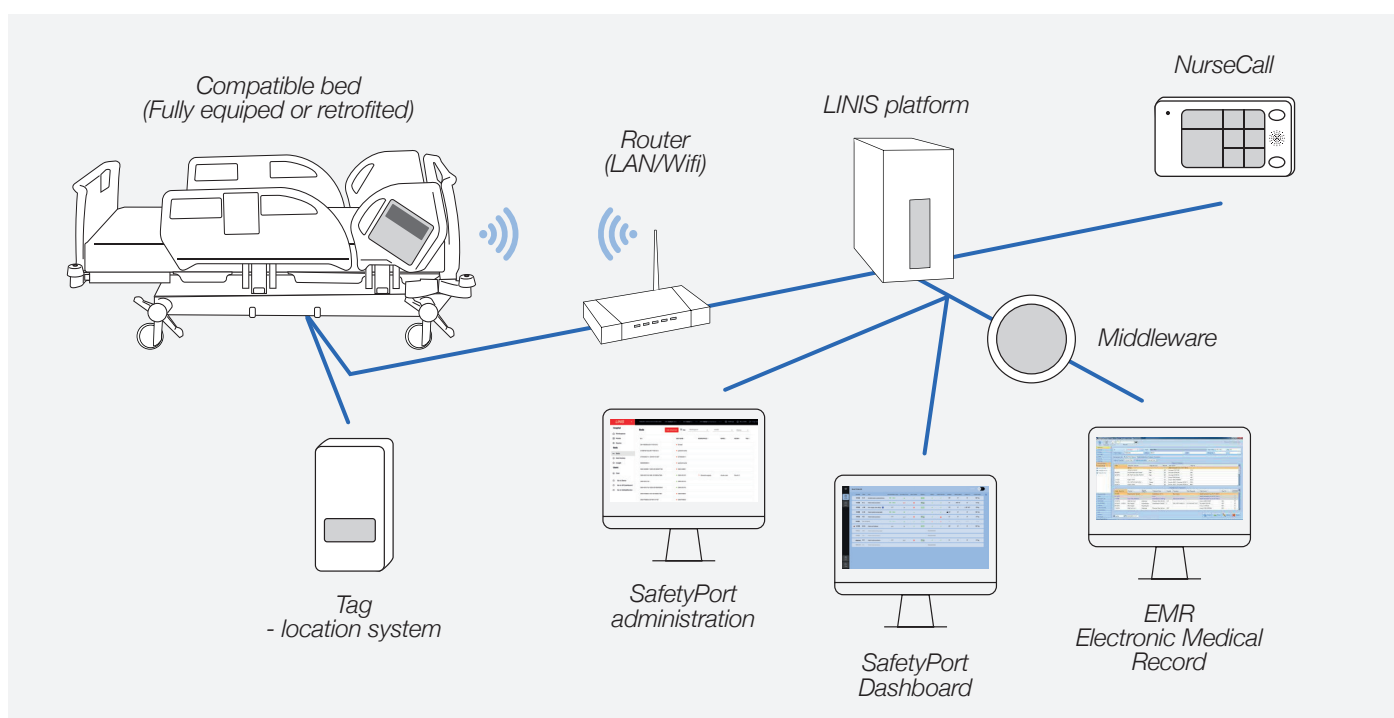

IT OPTIONS

Server operating system	Linux OS (preferred), Windows server (2016, 2019 and later versions)
Data format for export to the HIS	HL7 v2, xml, webservices
Graphical interface for display at the nurse station	Microsoft Edge, Mozilla Firefox, Google Chrome, Opera, Safari

wissner-
bosserhoff



Members of LINET Group

LINET spol. s r.o.

Želčevčice 5 | 274 01 Slaný | Czech Republic

tel.: +420 312 576 400 | fax: +420 312 522 668 | e-mail: info@linet.com | www.linet.com